



Paint 3D

Universal coloring system.  
3D effect pastes.

Universal coloring system.  
3D effect pastes.

# Paint 3D

**Paint 3D is a universal coloring system for the individual staining and glazing of dental ceramics in the CTE range from  $9.4$  to  $17.5 \times 10^{-6}/K$ .**

Ideally suited for milling, pressing and veneering ceramics – including zirconium oxide, lithium disilicate and other glass ceramics – Paint 3D offers a wide range of applications and maximum flexibility in everyday laboratory work.

The innovative 3D effect pastes open up new creative possibilities and give surfaces an impressive three-dimensional depth and texture. Whether for subtle accents or striking effects: These pastes make it possible to transform even the smallest details into realistic dental works of art.

Thanks to their special consistency, the pastes can be applied and modeled effortlessly, making them ideal for demanding projects where precision and efficiency are paramount. Paint 3D supports time-saving workflows without compromising on esthetics or attention to detail – for professional results at the highest level.



## Technical data



**Typ**  
Low-fusing glass  
ceramic



**Pastes**  
Self-glazing  
3D effect pastes



**CTE range**  
 $9.4 - 17.5 \times 10^{-6}/K$



**Applications**  
Zirconium oxide,  
Lithium-Disilicate,  
Veneering ceramics



**Application techniques**  
Wet-in-Wet,  
Classic

## Highlights


---

### 1 Universal coloring system

Versatile application on dental ceramics in the CTE range of  $9.4$  to  $17.5 \times 10^{-6}/K$



Applicable to zirconia, milled ceramics, pressed ceramics, and layered ceramics (full-ceramic & metal-ceramic)

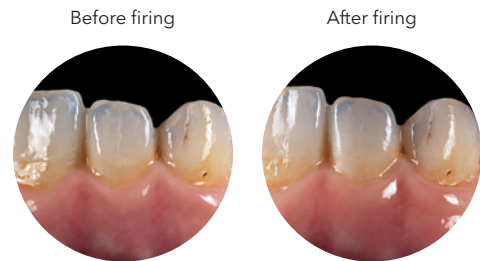

Firing-stable colors – what you see is what you get!

### 2 3D effect pastes

Creates three-dimensional effects and textures for a lively, natural esthetics



Unique optical depth, even at minimal layer thickness



Special texture paste for voluminous, three-dimensional surface textures

### 3 Easy handling

Intuitive and effortless application for time-saving and consistent results



Ready-to-use consistency



Self-glazing properties for easy application



# Delivery forms

Colors  
White Esthetics - pastes with fluorescence

Dentin



Paste per jar  
4g

Enamel



Effect



Intense Effect



Glaze & Texture



Pink Esthetics - pastes without fluorescence

Base



Effect



Glaze & Texture



Universal Liquid



767512

Starter Kit



767479

Gingiva Kit



767512

## Manufacturer

DSSM AG  
Im alten Riet 9  
9494 Schaan  
Liechtenstein  
sagemax.com